



Standard Features


- ▲ Checks bracelet for proper operation
- ▲ Remembers last bracelet test via recall button
- ▲ Identifies end-of-service bracelets
- ▲ Back-lit display
- ▲ Multi-color test result indicator
- ▲ Tests door coverage and door alarms
- ▲ Equipped with intelligent ambient light-sensing technology
- ▲ Built-in low battery warning
- ▲ Designed for a minimum of 2 years service

Specifications

- ▲ Dimensions: 4" x 2" x .75"
- ▲ Weight: 2.4 oz
- ▲ Power requirements:
Long life/non-replaceable battery

WanderGuard[®] Universal Tester Model 18029

The Universal Tester has been designed to provide a simple and effective way to test your WanderGuard ID System. It is a wireless portable device for testing WanderGuard ID signaling device bracelets and for testing working operation for WanderGuard ID door monitors and door alarms.

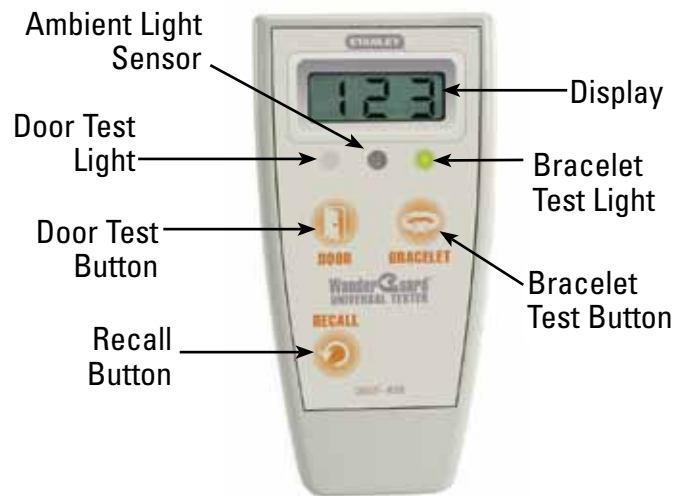

STANLEY[®]

Healthcare Solutions

Toll-Free: 800-824-4494

Office: 402-475-4002

Fax: 402-475-4281

www.stanleyhealthcare.com